

# AED Plus®

## Technical Specifications



### Measuring Chest Compression with Real CPR Help

The force required to depress the chest during CPR varies with the patient's size and build. The standard measure of chest compression quality, however, is not force but depth. The Real CPR Help® technology in ZOLL's CPR-D-padz® includes a hand-placement locator, an accelerometer, electronics, and a processing algorithm that work together to measure vertical displacement in space as each compression occurs.

### One Size Fits All


ZOLL's one-piece electrode design accounts for anatomical variation. Based on extensive human data, CPR-D-padz meet the anthropometric chest characteristics of 99% of the population. If needed, the lower (apex) electrode can be separated and adjusted to accommodate the remaining 1% of the population.

### Simplified Placement


Affixing two separate electrode pads to the patient's bare chest can be confusing to a lay rescuer. ZOLL's CPR-D-padz simplify this step by guiding placement of the red crosshairs at the center of the imaginary line connecting the patient's nipples. Once in place, the hand-locator and the two electrode pads fall naturally into optimal position for both defibrillation and CPR.

### Five-Year Shelf Life

All AED electrodes transmit defibrillating electricity into the patient via metal in close contact with a salt-infused gel that is positioned between the metal and the skin. Over time, however, the salt in the gel will corrode the metal and eventually compromise electrode functionality. ZOLL's novel electrode design includes a sacrificial element that prevents significant corrosion for five years, which is unmatched in the market today.

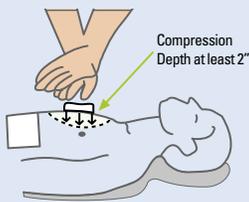

Real CPR Help® provides unique assistance to rescuers with real-time feedback on CPR compression depth and rate.

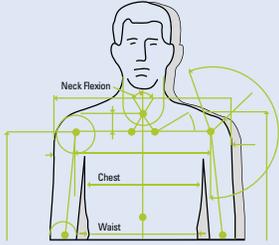

ZOLL's one-piece CPR-D-padz is designed to fit 99% of the population's chest anatomy.



CPR-D-padz offers clear anatomical placement illustrations and a CPR hand positioning landmark.



CPR-D-padz comes complete with rescue essentials, including a barrier mask, a razor, scissors, disposable gloves, and a towelette.

## AED Plus Specifications

### Defibrillator

**Waveform:** Rectilinear Biphasic

**Defibrillator Charge Hold Time:** 30 seconds

**Energy Selection:** Automatic preprogrammed selection (120J, 150J, 200J)

**Patient Safety:** All patient connections are electrically isolated

**Charge Time:** Less than 10 seconds with new batteries

**Electrodes:** ZOLL CPR-D-padz, pedi-padz® II or stat-padz® II

**Self-test:** Configurable automatic self-test from 1 to 7 days. Default = every 7 days. Monthly full-energy test (200J).

**Automatic Self-Test Checks:** Battery capacity, electrode connection, electrocardiogram and charge/discharge circuits, microprocessor hardware and software, CPR circuitry and CPR-D sensor, and audio circuitry

**CPR: Metronome Rate:** Variable 60 to 100 CPM

**Depth:** 3/4" to 3.5"; 1.9 to 8.9 cm

**Defibrillation Advisory:** Evaluates electrode connection and patient ECG to determine if defibrillation is required

**Shockable Rhythms:** Ventricular fibrillation with average amplitude >100 microvolts and wide complex ventricular tachycardia with rates greater than 150 BPM for adults, 200 BPM for pediatrics. For ECG Analysis Algorithm sensitivity and specificity, refer to AED Plus Administrator's Guide.

**Patient Impedance Measurement Range:** 0 to 300 ohms

**Defibrillator:** Protected ECG circuitry

**Display Format:** Optional LCD with moving bar

**Size:** 2.6" x 1.3"; 6.6 cm x 3.3 cm

**Optional ECG Viewing Window:** 2.6 seconds

**Optional ECG Display Sweep Speed:** 25 mm/sec; 1"/sec

**Battery Capacity:** Typical new (20°C) = 5 years (225 shocks) or 13 hours continuous monitoring. End of life designated by Red X (typical remaining shocks = 9).

**Data Recording and Storage:** 50 minutes of ECG and CPR data. If audio recording option is installed and enabled, 20 minutes of audio recording, ECG, and CPR data. If audio recording is disabled, 7 hours of ECG and CPR data.

**PC Minimum Requirements For Configuration and Patient Data Recovery:** Windows® 98, Windows® 2000, Windows® NT, Windows® XP, Windows® 7, IBM-compatible PII with 16550 UART (or higher) computer. 64MB RAM. VGA monitor or better. CD-ROM drive. IrDA port. 2MB disk space.

### Device

**Size:** (H x W x D) 5.25"x 9.50" x 11.50"; 13.3 cm x 24.1 cm x 29.2 cm

**Weight:** 6.7 lbs; 3.1 kg

**Power:** User-Replaceable Batteries. 10 Type 123A Photo Flash lithium manganese dioxide batteries

**Device Classification:** Class II and internally powered per EN60601-1

**Design Standards:** Meets applicable requirements of UL 2601, AAMI DF-39, IEC 601-2-4, EN60601-1, IEC60601-1-2

### Environmental

**Operating Temperature:** 32° to 122°F; 0° to 50°C

**Storage Temperature:** -22° to 140°F; -30° to 60°C

**Humidity:** 10 to 95% relative humidity, non-condensing

**Vibration:** MIL Std. 810F, Min. Helicopter Test

**Shock:** IEC 68-2-27; 100G

**Altitude:** -300 to 15,000 ft.; -91 m to 4573 m

**Particle and Water Ingress:** IP-55

### CPR-D-padz

**Shelf Life:** 5 years

**Conductive Gel:** Polymer Hydrogel

**Conductive Element:** Tin

**Packaging:** Multilayer foil laminate pouch

**Impedance Class:** Low

**Cable Length:** 48 in (1.2 m)

**Sternum:** Length: 6.1 in (15.5 cm); Width: 5.0 in (12.7 cm); Length, conductive gel: 3.5 in (8.9 cm); Width, conductive gel: 3.5 in (8.9 cm); Area, conductive gel: 12.3 sq in (79.0 sq cm)

**Apex:** Length: 6.1 in (15.5 cm); Width: 5.6 in (14.1 cm); Length, conductive gel: 3.5 in (8.9 cm); Width, conductive gel: 3.5 in (8.9 cm); Area, conductive gel: 12.3 sq in (79.0 sq cm)

**Complete Assembly:** Folded Length: 7.6 in (19.4 cm); Folded width: 7.0 in (17.8 cm); Folded height: 1.5 in (3.8 cm)

**Design Standards:** Meets applicable requirements of ANSI/AAMI/ISO DF-39-1993

**ADVANCING  
RESUSCITATION.  
TODAY.®**

ZOLL Medical Corporation  
Worldwide Headquarters  
269 Mill Road  
Chelmsford, MA 01824  
978-421-9655  
800-348-9011

For subsidiary addresses and fax numbers, as well as other global locations, please go to [www.zoll.com/contacts](http://www.zoll.com/contacts).

Specifications subject to change without notice.

©2011 ZOLL Medical Corporation. All rights reserved. Advancing Resuscitation Today, AED Plus, CPR-D-padz, pedi-padz II, Real CPR Help, stat-padz, and ZOLL are trademarks or registered trademarks of ZOLL Medical Corporation in the United States and/or other countries. All other trademarks are the property of their respective owners.

Printed in U.S.A. 121210 9656-0158

**ZOLL**®